

# Mapping of RNA modifications by direct Nanopore sequencing and JACUSA2

Christoph Dieterich<sup>\*1,2,3</sup>, Amina Lemsara<sup>1,2</sup>, and Isabel Naarmann-de Vries<sup>1,2,3</sup>

<sup>1</sup>Klaus Tschira Institute for Integrative Computational Cardiology, University Heidelberg, 69120 Heidelberg, Germany

<sup>2</sup>Department of Internal Medicine III (Cardiology, Angiology, and Pneumology), University Hospital Heidelberg, 69120 Heidelberg, Germany

<sup>3</sup>German Centre for Cardiovascular Research (DZHK)-Partner Site Heidelberg/Mannheim, 69120 Heidelberg, Germany

## Abstract

to be written

**Keywords:** Bayesian, 10X Genomics, Cell barcode assignment, Nonsense-mediated mRNA decay (NMD)

## INTRODUCTION

Chemical modifications on DNA and histones, also known as epigenetics marks, strongly impact gene expression during cell differentiation and in several other biological programs. In the 1970s, it was recognized that RNA is also subjected to extensive covalent modification, and studies in the late 1980s revealed the widespread deamination of bases (termed RNA editing), which can lead to recoding if it occurs within coding sequences. Impressive development in the RNA modification field occurred during the past eight years, with the discovery of an extensive layer of base modifications in mRNAs. These can influence gene expression and have been already shown to be involved in primary cellular programs such as stem cell differentiation, response to stress, and the circadian clock. The study of RNA modifications and their effects is now referred to as epitranscriptomics, and it reveals striking similarities to what is known for epigenomics. To date thirteen distinct modifications have been identified on mRNA transcripts [?]. These modifications are catalyzed by a variety of dedicated enzymes and can be divided into two classes: modifications of cap-adjacent nucleotides and internal modifications.

---

\*christoph.dieterich@uni-heidelberg.de

32 In contrast to the m7G cap, the impact of internal modifications on gene  
 33 regulation has been less studied apart from RNA editing, which is mediated  
 34 by RNA deaminases (e.g. the ADAR family). The most widespread internal  
 35 mRNA modification is N6-methyladenosine (m6A). By modulating the pro-  
 36 cessing of mRNA, m6A can regulate a wide range of physiological processes  
 37 and its alteration has been linked to several diseases [???]. The modifica-  
 38 tion is catalyzed co-transcriptionally by a Mega-Dalton methyltransferase  
 39 complex, which includes the heterodimer METTL3-METTL14 and other  
 40 associated subunits [?]. This modification is reversible since two proteins of  
 41 the AlkB-family demethylases can remove m6A from mRNA transcripts [??].  
 42 In mammals, m6A preferentially localizes within long internal exons and at  
 43 the beginning of terminal exons at so-called DRACH motif (D = A/G/U,  
 44 R = A/G, H = A/C/U) sites [???]. Once deposited, m6A is recognized  
 45 by several reader proteins that can affect the fate of mRNA transcripts in  
 46 nearly every step of the mRNA life cycle, which includes alternative splicing  
 47 [?]. The best-described readers are the YTH domain family of proteins that  
 48 decode the signal and mediate m6A functions. By affecting RNA structure,  
 49 m6A can also indirectly influence the association of additional RNA-binding  
 50 proteins (RBPs) and the assembly of larger messenger ribonucleoprotein  
 51 particles (mRNPs).

52 Several approaches have been presented to map RNA modifications on  
 53 RNA. Herein, we focus on mRNA modification site detection in general and  
 54 on m6A in particular where antibody-based protocols (miCLIP), methylation-  
 55 sensitive restriction enzyme assays (MazF) or transgenic approaches (TRIBE,  
 56 DART) have been presented. All of the aforementioned approaches rely on  
 57 high-throughput sequencing on the Illumina platform. This typically in-  
 58 volves cDNA synthesis by reverse transcription and PCR-based library am-  
 59 plification. One recent addition to the tool is direct RNA single molecule  
 60 sequencing on the Oxford Nanopore Technology platform. While our software  
 61 workflow is able to deal with Illumina and Nanopore-based approaches, the  
 62 latter is the principal topic of our methods article.

## 63 MATERIALS

### 64 ONT direct RNA sequencing

- 65 1. 500 ng polyA<sup>+</sup> RNA isolated from total RNA e.g. with Oligotex  
 66 mRNA kit (Qiagen) or Dynabeads oligo dT<sub>25</sub> beads (Thermo Fisher  
 67 Scientific). Store RNA at -80 °C and the mRNA purification kit as  
 68 recommended by the manufacturer.
- 69 2. Nuclease-free water. Store at room temperature.
- 70 3. Direct RNA-sequencing kit (SQK-RNA002, Oxford Nanopore Tech-  
 71 nologies). Store at -20 °C.

- 72 4. NEBNext Quick Ligation Reaction Buffer (New England Biolabs).  
73 Store at -20 °C.
- 74 5. T4 DNA Ligase (New England Biolabs). Store at -20 °C.
- 75 6. dNTP Mix (10 mM each). Store at -20 °C.
- 76 7. SuperScript IV Reverse Transcriptase (Thermo Fisher Scientific). Store  
77 at -20 °C.
- 78 8. Agencourt RNAClean XP beads (Beckman Coulter). Store at 4 °C.
- 79 9. 70 % ethanol, freshly prepared.
- 80 10. Qubit dsDNA HS assay kit and Qubit Fluorometer (Thermo Fisher  
81 Scientific).
- 82 11. Flow cell priming kit (EXP-FLP002, Oxford Nanopore Technologies).  
83 Store at -20 °C.
- 84 12. Thermo cycler.
- 85 13. Gentle rotator mixer.
- 86 14. Magnetic stand for 1.5 ml tubes.
- 87 15. 1.5 ml DNA LoBind tubes (Eppendorf), 0.2 ml PCR tubes.
- 88 16. MinION or GridION sequencing device and MinION R9.4.1 Flow cells  
89 (FLO-MIN106D, Oxford Nanopore Technologies). Store Flow cells at  
90 4 °C.

## 91 Preparation of an *in vitro* transcriptome sample

92 Isabel

## 93 Hardware requirements

94 All analyses have been performed/tested on two alternative hardware sys-  
95 tems: a standard Linux desktop computer or an Apple iMac (Retina 5K,  
96 ultimo 2014). The workflow requires a multi-core processor system with  
97 minimal main memory of 16GB RAM and several GBs of free disk space  
98 (depending on data set size).

## 99 Software dependencies and installation

100 Our analysis workflow has few requirements, which are detailed in Table 2.  
101 Specifically, to execute our workflow, the following prerequisites are neces-  
102 sary: a BASH shell, a JAVA runtime environment, a working PERL and  
103 R installation. Additional i.e. non-standard software to process and map  
104 Nanopore reads (bedtools, samtools and Minimap2) are obligatory, while  
105 the installation of a Nanopore read simulator (NanoSim) is optional and de-  
106 pends on your use case. Table ?? lists some additional R packages, which are  
107 required to run the R code. Detailed instructions on how to setup are found  
108 under [https://github.com/dieterich-lab/MiMB\\_JACUSA2\\_chapter](https://github.com/dieterich-lab/MiMB_JACUSA2_chapter)

## 109 METHODS

110 Overview Figure 1

### 111 Nanopore direct RNA sequencing

112 Isabel

- 113 1. Adjust 500 ng polyA<sup>+</sup> RNA to a total volume of 9  $\mu$ l with nuclease-  
114 free water. Complete RT adapter ligation reaction (in 0.2 ml PCR  
115 tube) with 3  $\mu$ l NEBNext Quick Ligation Reaction Buffer, 0.5  $\mu$ l  
116 RNA CS (RCS, from SQK-RNA002), 1  $\mu$ l RT-Adapter (RTA, from  
117 SQK-RNA002) and 1.5  $\mu$ l T4 DNA Ligase. Incubate 10 min at room  
118 temperature.
- 119 2. Prepare reverse transcription master mix during ligation: 9  $\mu$ l nuclease-  
120 free water, 2  $\mu$ l 10 mM dNTPs, 8  $\mu$ l 5x SuperScript IV first strand  
121 buffer, 4  $\mu$ l 0.1 mM DTT.
- 122 3. Add the reverse transcription master mix to the ligation reaction and  
123 mix by pipetting. Add 2  $\mu$ l SuperScript IV reverse transcriptase and  
124 mix by pipetting. Incubate in a thermocycler with the following pro-  
125 tocol: 50 min at 50 °C, 10 min at 70 °C, cool down to 4 °C.
- 126 4. Let the Agencourt RNAClean XP beads come to room temperature  
127 during incubation of reverse transcription. Carefully resuspend beads  
128 before use. Transfer reaction to a 1.5 ml DNA LoBind tube and mix  
129 with 72  $\mu$ l Agencourt RNAClean XP beads. Incubate 5 min at room  
130 temperature on a gentle rotator mixer.
- 131 5. Collect beads on a magnetic stand and remove supernatant. Wash  
132 pelleted beads two times (30 sec) with 200  $\mu$ l freshly prepared 70 %  
133 ethanol. Remove supernatant. Spin sample down and place on magnet  
134 again. Remove any residual ethanol.

- 135 6. Resuspend beads in 20  $\mu$ l nuclease-free water by gentle flicking and  
136 incubate 5 min at room temperature on a gentle rotator mixer. Collect  
137 beads on a magnetic stand and transfer 20  $\mu$ l eluate in a fresh 1.5 ml  
138 DNA LoBind tube.
- 139 7. For ligation of the RMX adapter, add the following to 20  $\mu$ l eluate: 8  
140  $\mu$ l NEBNext Quick Ligation Reaction Buffer, 6  $\mu$ l RMX (from SQK-  
141 RNA002), 3  $\mu$ l nuclease-free water, 3  $\mu$ l T4 DNA Ligase. Mix by  
142 pipetting and incubate 10 min at room temperature.
- 143 8. Add 40  $\mu$ l carefully resuspended Agencourt RNAClean XP beads to  
144 the reaction and mix by pipetting. Incubate 5 min at room tempera-  
145 ture on a gentle rotator mixer.
- 146 9. Collect beads on a magnetic stand and remove supernatant. Wash  
147 pelleted beads two times with 150  $\mu$ l wash buffer (WSB, from SQK-  
148 RNA002). Resuspend beads by flicking, spin down and return to mag-  
149 netic stand. Remove supernatant from pelleted beads.
- 150 10. Resuspend beads in 21  $\mu$ l elution buffer (EB, from SQK-RNA002) by  
151 gentle flicking and incubate 5 min at room temperature on a gentle  
152 rotator mixer. Pellet beads on a magnetic stand and transfer 21  $\mu$ l  
153 eluate in a fresh 1.5 ml DNA LoBind tube.
- 154 11. Quantify 1  $\mu$ l of the library on a Qubit fluorometer with the Qubit  
155 dsDNA HS kit according to the manufacturerers protocol. Concentra-  
156 tion should be usually in the range of 5 - 10 ng/ $\mu$ l.
- 157 12. Insert MinION R9.4.1 Flow cell in the MinION or GridION sequenc-  
158 ing device and perform Flow cell check in the MinKNOW software.  
159 For successful sequencing of mammalian polyA<sup>+</sup> RNA at least 1,000  
160 available pores are recommended.
- 161 13. Prepare Priming Mix by adding 30  $\mu$ l flush tether (FLT, from EXP-  
162 FLP002) to a vial of flush buffer (FB, from EXP-FLP002) and mix by  
163 pipetting. Open priming port. Remove air bubble from priming port  
164 by inserting the tip of a P1000 pipette into the priming port and slowly  
165 dialing up, until a small volume of storage buffer enters the pipette  
166 tip. Load 800  $\mu$ l Priming Mix via the priming port and carefully avoid  
167 introduction of air bubbles. Close the priming port and wait for 5 min.
- 168 14. Mix 20  $\mu$ l library with 17.5  $\mu$ l nuclease-free water and 37.5  $\mu$ l RNA run-  
169 ning buffer (RRB, from SQK-RNA002) and mix by pipetting. Open  
170 the priming port and the sample port. Load 200  $\mu$ l Priming Mix via  
171 the priming port. Mix library by pipetting just before loading and  
172 load dropwise via the sample port. Carefully avoid introduction of air  
173 bubbles. Close the sample port and the priming port.

174 15. Start sequencing for 48 to 72 h in the MinKNOW software. Choose  
175 direct RNA-sequencing kit and high-accuracy basecalling as param-  
176 eters. We recommend to adjust the output filter to a minimum Q score  
177 of 7 (instead of 9).

## 178 Preparation of an *in vitro* transcriptome sample

179 Isabel

## 180 Nanopore read processing

181 Minimap2 and samtools

182 Christoph

## 183 Use Case 1: Comparison of wildtype and knock-out samples

184 Xpore data

185 Christoph

## 186 Use Case 2: Comparison of wildtype and IVT samples

187 Christoph

## 188 Use Case 3: Comparison of wildtype to simulated IVT sample

189 Christoph

## 190 NOTES

### 191 Tips and Tricks

## 192 ACKNOWLEDGMENTS

193 The authors would like to thank Etienne Boileau, Thiago Britto Borges,  
194 Tobias Jakobi for proof-reading and comments. The authors are grateful  
195 to Marek Franitza for running the experiments on the 10x platform and to  
196 Christian Becker for running ONT sequencing. This work was supported by  
197 Informatics for Life funded by the Klaus Tschira Foundation.

## FIGURE CAPTIONS

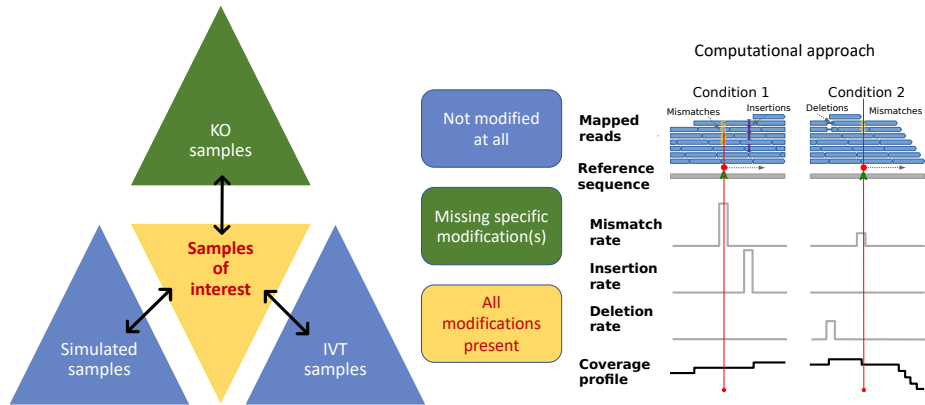


Figure 1: **General outline of RNA modification detection by JACUSA2.** A key feature of our approach is that multiple replicates can be compared as shown on the left. Samples of interests where all modifications are present could be compared with either KO samples where the modification of interest is missing or IVT/simulated samples where all modifications are absent. Read stacks (in blue) are compared head-to-head as shown on the right.

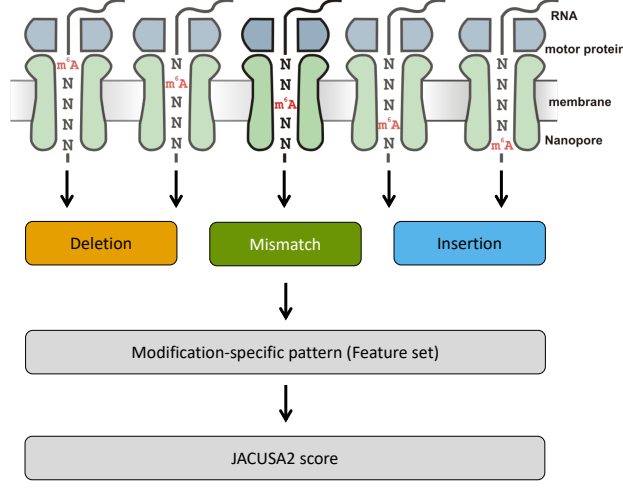


Figure 2: **Motivation of 5mer context for RNA modification mapping.** The nanopore covers 5 consecutive RNA residues. That is why we consider a 5mer context and derive 3 principal features for every position within a given 5 mer (15 features in total, with a central A residue in this example). We evaluate each feature set by previously learned patterns and compute a final score for modification site detection.



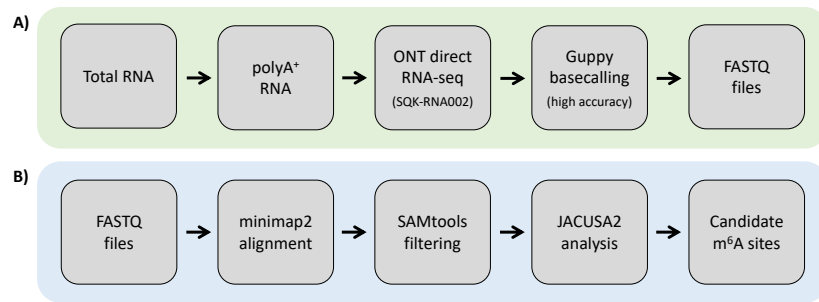


Figure 3: **Experimental and computational workflow.** tbd

Software	Version	Description
Minimap2	<a href="https://github.com/lh3/minimap2">https://github.com/lh3/minimap2</a> v2.22 or later	<a href="https://lh3.github.io/minimap2/">https://lh3.github.io/minimap2/</a>
samtools	<a href="https://github.com/samtools/samtools">https://github.com/samtools/samtools</a> v1.12 or later	<a href="http://samtools.github.io/">http://samtools.github.io/</a>
JAVA	openjdk 11.0.12 2021-07-20 - JAVA 11 or later	OpenJDK Runtime Environment
R	<a href="https://www.r-project.org/">https://www.r-project.org/</a> version 3.5.1 or later	The R Project for Statistical Computing
PERL	<a href="https://www.perl.org/">https://www.perl.org/</a> version 5.28.1 or later	Perl is a highly capable, feature-rich programming language
BASH, sed, awk	should be part of your Linux distribution	Misc.
bedtools	<a href="https://github.com/arq5x/bedtools2">https://github.com/arq5x/bedtools2</a> version 2.29.2 or later	Perl is a highly capable, feature-rich programming language
NanoSim	<a href="https://github.com/bcgsc/NanoSim">https://github.com/bcgsc/NanoSim</a> version 3.0.2 or later (optional)	NanoSim is a fast and scalable read simulator that captures the technology-specific features of ONT data

Table 1: **Software dependencies** blubba

## 199 TABLE CAPTIONS

## 200 TABLES

R Pack- ages	Version	Description
ggplot2	<a href="https://cran.r-project.org/web/packages/ggplot2/index.html">https://cran.r-project.org/web/packages/ggplot2/index.html</a> - ggplot2_3.3.0 or later	ggplot2 is a system for declaratively creating graphics, based on The Grammar of Graphics.
NMF	<a href="https://cran.r-project.org/web/packages/NMF/index.html">https://cran.r-project.org/web/packages/NMF/index.html</a> - NMF_0.22.0 or later	Provides a framework to perform Non-negative Matrix Factorization (NMF).

Table 2: **R Package dependencies** blubba